

SOPRANO TITANIUM VERSION

Diode LASER HAIR REMOVAL

Even faster, Safty and Effective



14+ YEARS EXPERIENCE

AESTHETICS MEDICAL - EQUIPMENT

OEM & ODM



100+

Patented technologies



500+

OEM real case



45+

R&D team



10000+

World-class supplier resources



Introducing Soprano Titanium A breakthrough hair removal platform that combines ultimate functionality with unparalleled comfort. Thanks to its 2 different spot size and advanced cooling system, treatments are now much faster and pain free- offering a significant improved solution for patients and practitioners alike. Soprano Titanium leverages the benefits of three combined wavelengths, with an improved patient experience and a business-oriented approach, creating a unique and result-driven new solution in the world of professional hair removal.

The Titanium's System Highlights

✓ Spot Size for Optional

15*40mm, 12*30mm

With 2 different spot sizes, it is able to work on large areas as well as small areas. This breakthrough applicator makes treatment more convenient and virtually painless, significantly reducing treatment time and ensuring a more comfortable experience for patients.

✓ USA Coherent Laser

1600W or 1200W

All kinds of power and wavelength laser module can reach your requirements



✓ ICE Plus Technology

Continuous Contact Cooling

Continuous Contact Cooling ICE Plus is an advanced technology that continuously cools skin and facilitates range-controlled cools skin and facilitates range-controlled temperature throughout the entire treatment,

The ICE Plus cold sapphire tip minimize esepidermal risks while maintaining heat within the dermis where hair follicles are treated, Temperature remains controlled and constant which contributes to the diode's efficiency. Patients can now enjoy a completely pain-free, cool and comfortable experience throughout the session.



✓ 15' Android Screen

3 intelligent Screen linkage

Friendly and Functional Equipped with a new large sized 15' Android screen, with pre-set parameters, practitioners can enjoy a more intuitive, user friendly work experience.




Synergetic Technologies. Superior Results.

SHR™ - SHR hair removal is the only clinically proven method of laser hair removal that is virtually painless, and it is the safest laser hair removal method for all skin types, including dark and tanned skin. Its unique gradual heating and rapid/fast pulse delivery method effectively damages the hair follicle and hinders re-growth, while preventing injury to the surrounding tissue.

SHR affects the melanin of the skin using a unique method, allowing the procedure to be performed on all skin types anytime of the year as well as on tanned skin.

3D technology

Featuring its exclusive 3D technology, Soprano Titanium combines the three most effective laser wavelengths into a single applicator, simultaneously targeting different tissue depths and anatomical structures within the hair follicle. By combining the absorption and penetration levels of different wavelengths, along with extended treatment coverage, increased comfort and low maintenance requirements, Soprano Titanium is the safest and most efficient hair removal treatment available today. 3D technology for one superior hair removal solution.

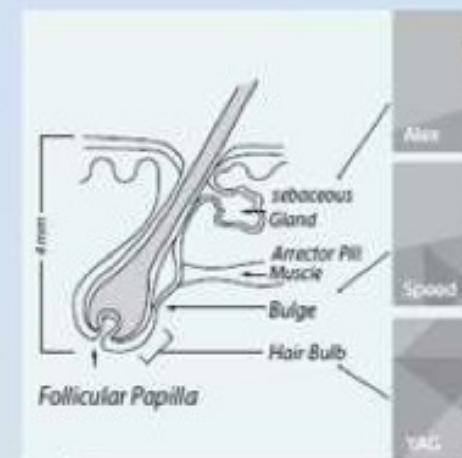

755nm

The 755nm wavelength, based on the Alex laser, offers a strong absorption by the chromophore melanin. This makes the 755nm an ideal wavelength in the removal of fine and light colored hair.

810nm

The 810nm wavelength is the gold standard for laser hair removal. It balances safety with efficacy by possessing good absorption in melanin and relatively deep penetration.

1064nm

The 1064nm wavelength is characterized by lower melanin absorption, making it a perfect solution for darker skin types. The 1064nm also offers deep penetration targeting the bulb and papilla of the hair follicle. With water and hemoglobin absorption, the 1064nm wavelength increases the thermal profile of the overall treatment, leading to highly effective hair removal, especially for deep hair follicles.

Three simultaneous wavelengths in one applicator



755nm

- Wavelength -

+

810nm

- Wavelength -

+

1064nm

- Wavelength -

Ideal for range of
hair colors

Safe for darker
skin types

Perfect solution for different
hair types and densities



Platform Benefits

- ⦿ **Golden laser bars from US Coherent with 30 million shots**
 - ⦿ **LARGE SPOT SIZE**
Different spot size for optional
 - ⦿ **COOLING SYSTEM**
Super sapphire touch directly cooling system ensures painless and comfortable treatment
 - ⦿ **TWO CONNECTORS**
 - ⦿ **CLINICALLY PROVEN 3D TECHNOLOGY**
 - ⦿ **INTUITIVE, USER FRIENDLY LARGE SCREEN**
(15'ANDROID) WITH PRE-SET PARAMETERS
 - ⦿ **I-VI SKIN TONES AND TANNED SKIN**
 - ⦿ **NO CONSUMABLES**
- 

Technical Parameters

Laser type	3 waves Diode Laser Hair Removal Soprano Titanium
Laser wavelength	755/810/1064nm+2nm
Spot size	15*40 mm ² , 13*20 mm ² (optional)
Pulse width	1-400 ms (adjustable) or 1-300ms
Energy	1-160J (adjustable), 280J/320J/380J optional
Frequency	1-10 Hz (adjustable)
Language	English, Spanish, Italian, German, French, Turkish or any other language as requested
Skin type	I-VI skin type
Rent System	Optional
Display	15' color 4K touch LCD display
Laser module power	600W/800W/1000W/1200W/1600W
Cooling	Water + Air + Semiconductor + A/C+TEC
Voltage	110V/220V+20V,50/60Hz

ZENLASE

Aesthetic Laser Equipment · Factory Direct

WHATSAPP

+86 151 0211 7919

EMAIL

info@zenlase.com

WEB

www.zenlase.com

HOURS

Mon-Sat 09:00-18:00 (GMT+8)

Quotations, orders, spare parts and after-sales for overseas buyers
are handled by the ZENLASE team.