

360° ANGLE SURROUNDING CRYOLIPOLYSIS SLIMMING MACHINE

3 CRYO+40K+BODY RF+FACE RF+LIPO LASER



7 in 1
MULTIFUNCTIONAL
COMBINATION



10.4" LARGE SCREEN



360° Refrigeration



Green light flashing



Exchangeable cryo
handles



Different sizes for
changeable



Adjustable time for
cold & hot

What is 360 cryolipolysis?


360-degree full range of cooling and heating, freezing at minus 10°C to positive 45°C heating, 4 groups cycle modes for operation, parameters can be set flexibly; 8 pcs (Changeable Cryo handle cups) different sizes cryo handles are suitable for different size areas and body shape.

NEW TECHNOLOGY---CRYO ANGLE

Cryo Angle is the latest fat freezing cooling technology which uses a special 360° applicator to target stubborn fat that is resistant to changes in diet and exercise, effectively freezing, destroying, and permanently eliminating the fat cells beneath the skin without damaging the surrounding layers.




New system

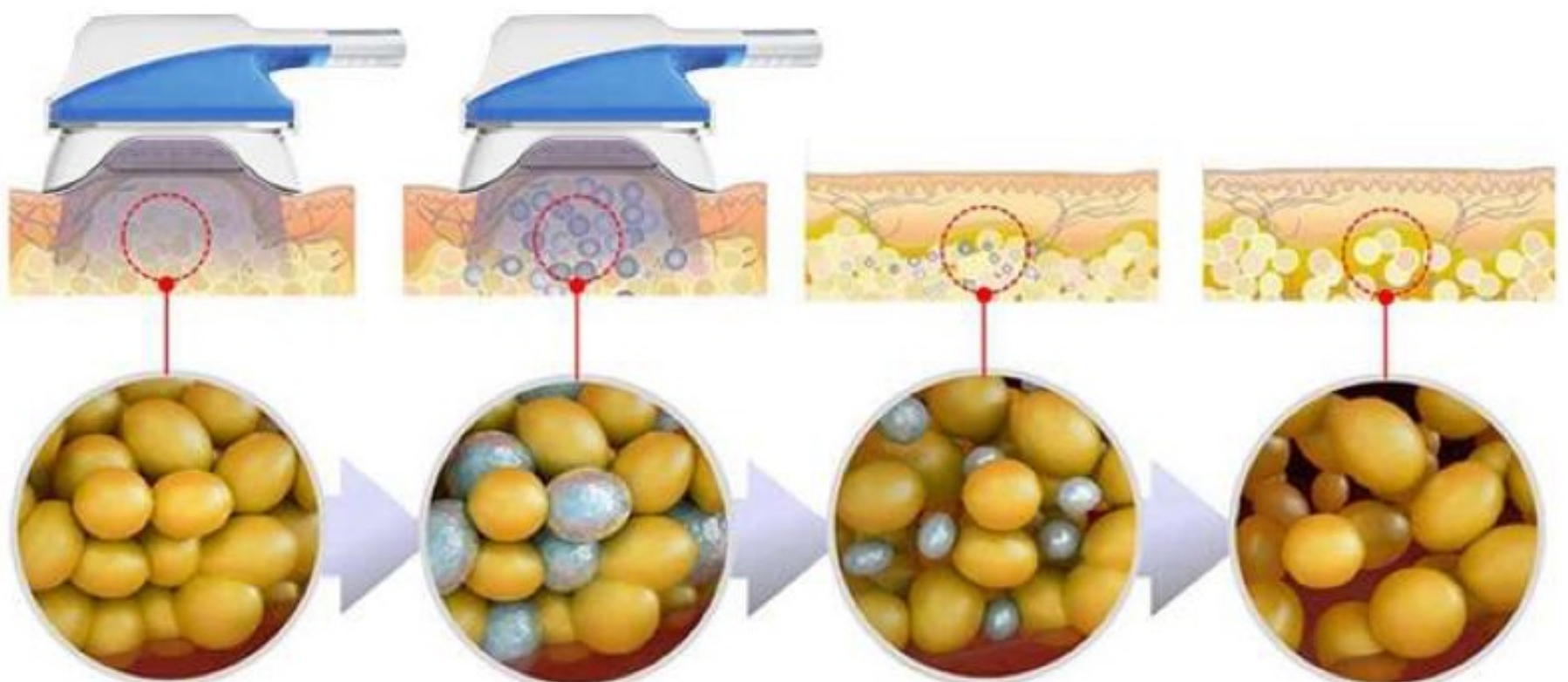

- Apply Coolshaping applicator at the target area, Start cryo-cooling to reach a target temperature.
- The fat cells in the target area begin responding to cryo-crystallization.
- Following the treatment, fat cells reach an biological death and removed by metabolism process.
- The fat cells are eliminated.

EXCHANGEABLE CRYO HANDLE CUPS

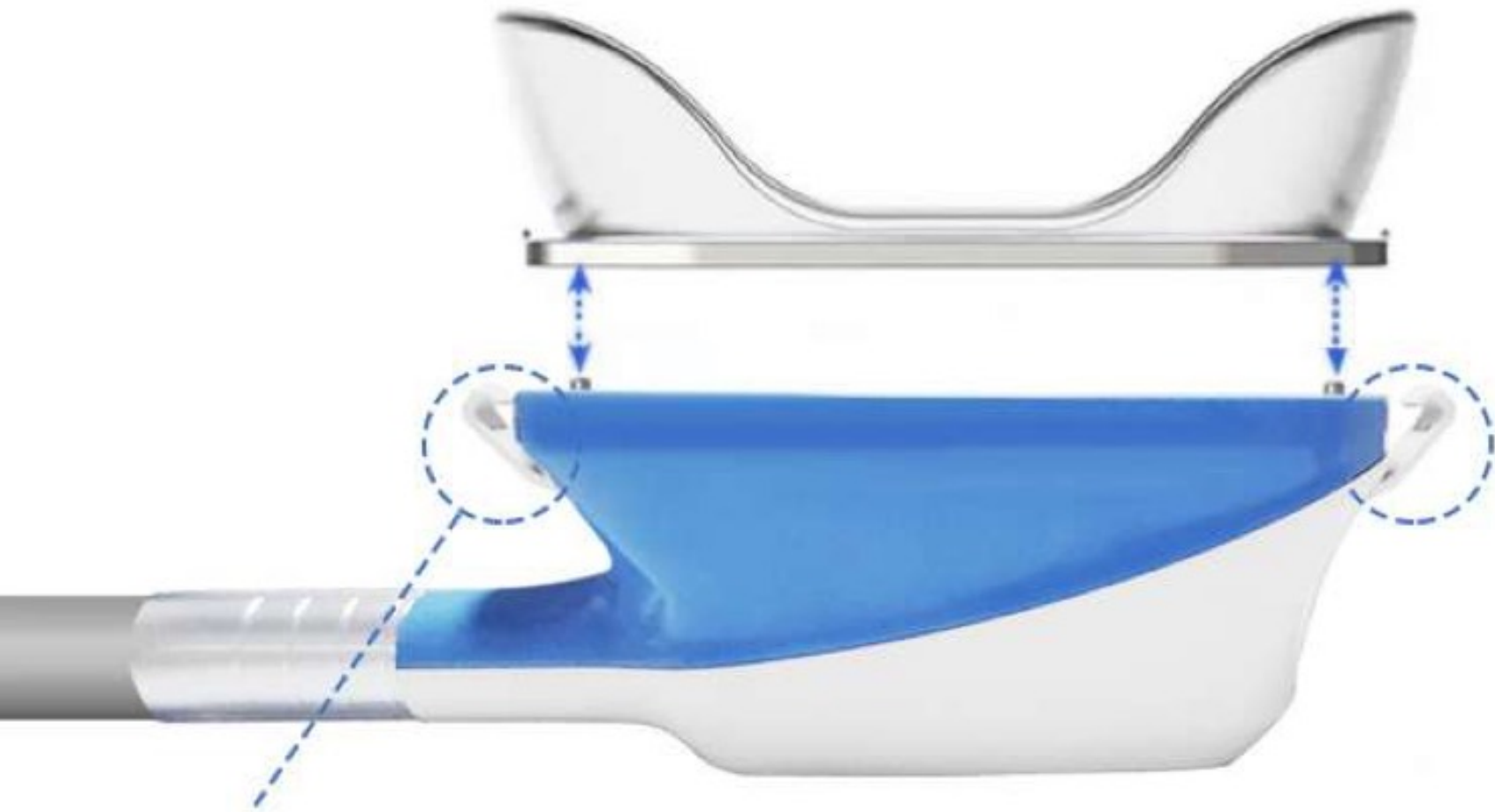

EASY TO CHANGE CRYO HANDLE CUPS

- CONTOURS ARE SECURELY FASTENED TO THE APPLICATOR
- QUICKLY CHANGE CONTOURS FOR YOUR NEXT TREATMENT

LARGE SIZE

Large cryo handle size: 1#195*85mm;
2#210*90mm;3#215*100mm;4#230*105mm

#3/#4 APPLICATION

for flanks.curved cup design allows for better placement and fit on curved or narrow parts of the body.

#2 APPLICATION

best matches with the contour of abdomen.most commonly used applicator used for this area.



MEDIUM SIZE


Small cryo handle size: 1#140*70mm;
2#155*80mm;3#165*90mm

#1 APPLICATION

ideal for longer vertical areas of fat(I.E.abdomen and flanks).also suggested for the arms.



COOLING CUP

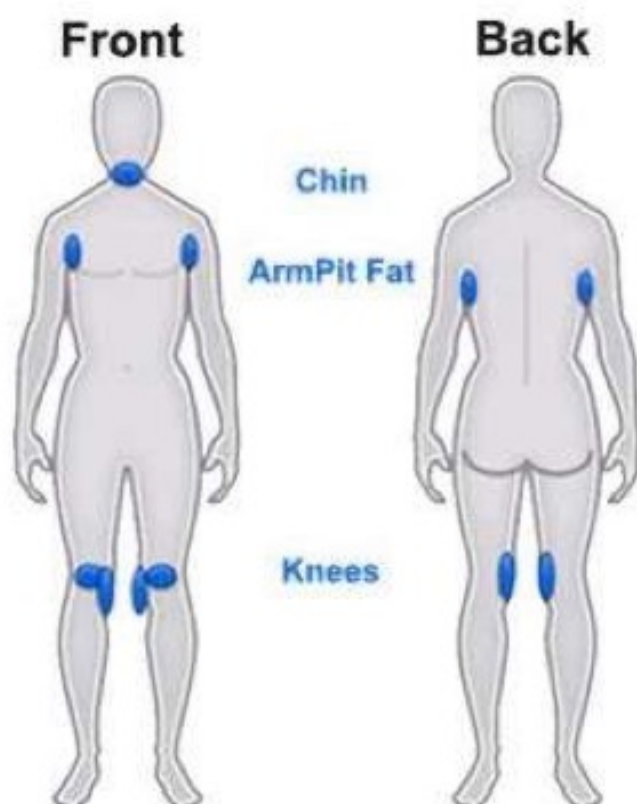

HANDLE SIZE: 80*40mm

Mini freeze

- ✦ Ideal for smaller bulges
- ✦ Safely targets delicate, hard-to-reach areas
- ✦ 35 minute cycle , Complete results in 1-2 visits
- ✦ A smooth, contoured, youthful jawline and chin
- ✦ No scalpels, injections, incisions, traumatic suction, or anesthesia



FreezeCurve

- Applicator of choice for flanks
- Curved cup design allows for better placement and fit on curved or narrow parts of the body



FreezeCore

- Best matches with the contour of abdomen
- Most commonly used applicator used for this area



FreezeFit

- Ideal for longer vertical areas of fat (i.e. abdomen and flanks)
- Also suggested for the arms and inner thighs (OUS)



Large size



Medium size



40K Cavitation

direct into the fatty layer, speedily vibrate deep-seated cellulite, produce numberless vacuum cavitation, mightily strike the fatty cells, let them produce inner cracking, and dissolve to be the free fatty acid.



Body RF

uses proprietary technology to create a highly efficient and tightly woven energy matrix. This dense energy matrix penetrates multiple layers of the dermis heating it from the inside out. Complete skin penetration ensures exponentially improved clinical efficacy and products clearly visible results.



Face RF

RF radio frequency is one kind of high frequency electromagnetic wave. It can penetrate the dermis and subcutaneous tissue. Face RF can be used for face lifting, wrinkle removal.



Lipo laser

utilizes the latest low level/cold laser technology to reduce girth and spot fat over almost all areas of the body without surgery, downtime or redness. A typical 40-minute waistline treatment can reduce an individual by $\frac{1}{2}$ to $\frac{3}{4}$ of an inch while a full protocol of eight treatments usually results in multiple inch loss.




INSTRUMENT INFORMATION

Model	QM-3S	Rated voltage	AC 110~220V
Function	7 in 1 CRYO	Rated power	1000W
Gross weight	39KG	Packaging	Wooden box
Frequency	50~60Hz	Package Size /L*H*W	64*47*64cm
Relative Humidity	≤ 80%		

TECHNICAL PARAMETERS

Handle control method	Touch	Large cryo handle size: 4#230*105mm	
Screen size	10.4	1#195*85mm 2#210*90mm 3#215*100mm	
Handle screen size	4.3	Small cryo handle size: 3#165*90mm	
Handle material	Silicon	1#140*70mm 2#155*80mm	
Handle temperature	-10°C/+45°C	Double chin handle size	80*40mm
Pressure output	0.08MPa	Cavitation head frequency	40KHz
Refrigeration method	Water cycle	Laser power	100mw/pcs
RF face	3MHz	Laser wavelength	650nm
RF frequency	5MHz	Number of laser lights	72PCS

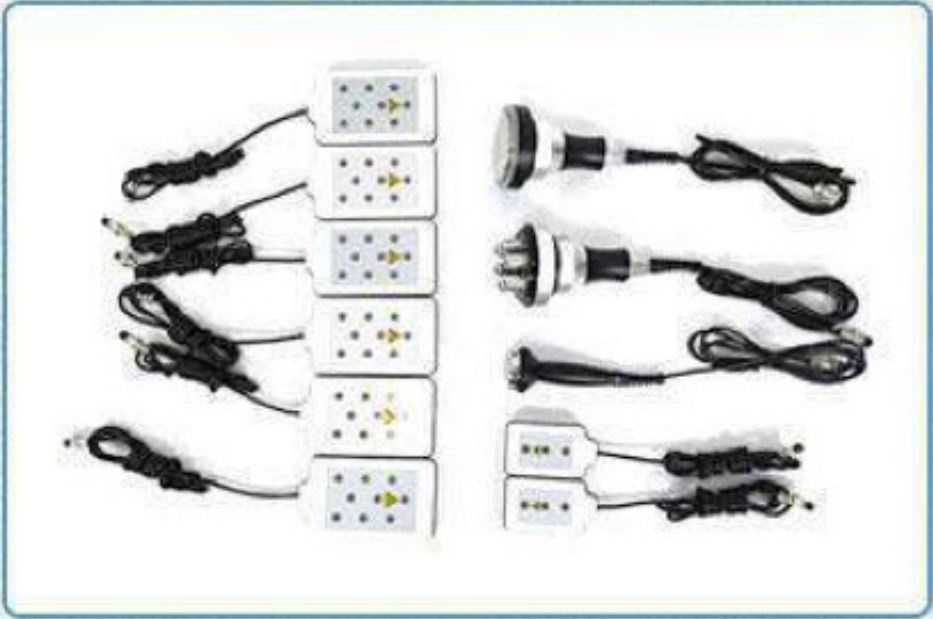

ZENLASE

Aesthetic Laser Equipment · Factory Direct

WHATSAPP

+86 151 0211 7919

EMAIL

info@zenlase.com

WEB

www.zenlase.com

HOURS

Mon-Sat 09:00-18:00 (GMT+8)

Quotations, orders, spare parts and after-sales for overseas buyers are handled by the ZENLASE team.