

New Air Pressure Shock Wave

QM-SW9+



What is Pneumatic Ballistic Shock Wave?

Pneumatic ballistic shock wave instrument is compressor generated pulse sound waves into precise ballistic, compare with the electromagnetic shock wave, more energy generated max to 6 bar. And also impact faster, configuration higher and more expensive.



Technical Parameters

Model	QM-SW9+
Frequency	1 - 21Hz
Energy	1 - 6Bar
Packing Size	58*46*38cm
Volume Weight	28kg
Power supply	AC100V-110V/AC220V-230V,50Hz/60Hz
Power	200W
Screen	8-inch touch screen
Treated Heads	12pcs including radial form, focus form and flat form

Advantages

Multiple waveforms
flat form/focus
form/radial form

Total 11pcs heads
more powerful
deeper penetration
Long lifetime

Easy to change
bullet

Max to 6 bar
max to 21Hz


Upgrade software
with ED therapy

Smarter programing
C and P mode
6 different preload


Easier
maintaining



The Interface

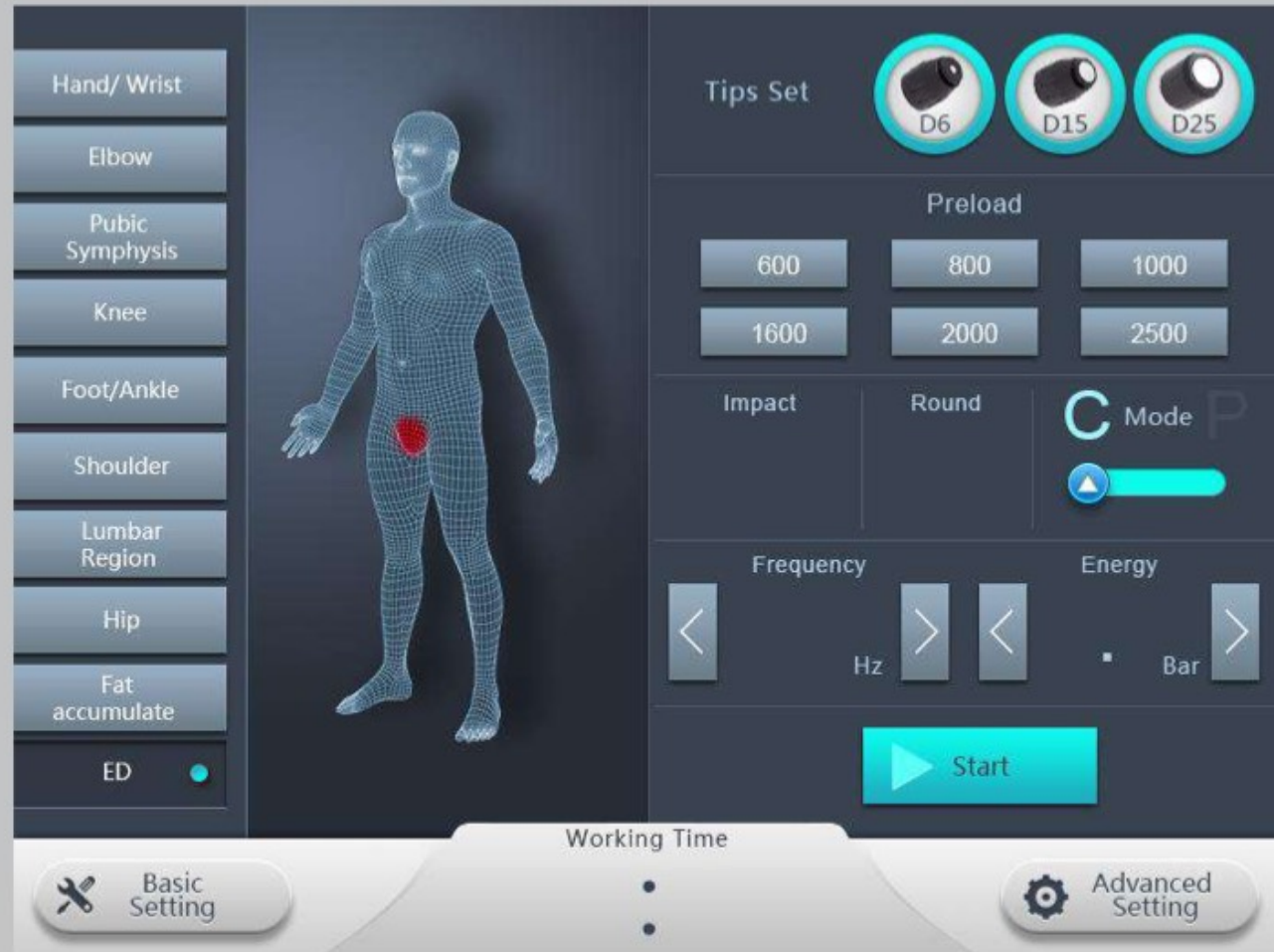

Handle and Tips

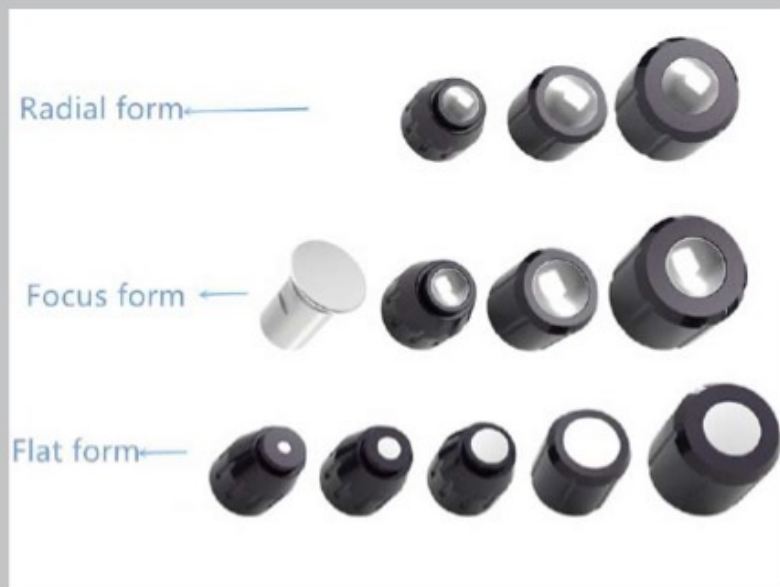

The upgraded handle and 12pcs tips
More impacting head choices
6- 8-10-15-20-25mm diameter heads

APPLICATIONS

01 Sports Injury Recovery

02 Rehabilitation Therapy

03 Physiotherapy Pain Relief

04 Cellulite Reduction

05 ED Treatment

Shock Wave Physiotherapy

Shockwave therapy is a multidisciplinary device used in orthopaedics, physiotherapy, sports medicine, urology and veterinary medicine.

Its main assets are fast pain relief and mobility restoration. Together with being a non-surgical therapy with no need for painkillers makes it an ideal therapy to speed up recovery and cure various indications causing acute or chronic pain.



Shock Wave Cellulite Treatment




Shock wave is a powerful and effective treatment in reducing cellulite, tightening and firming as well as smoothing of the skin.

shock wave can be used as a standalone application or in conjunction with other fat removal or skin tightening treatments to accelerate results.

Shock Wave for ED Therapy


Most men suffering from erectile dysfunction have vascular problems affecting the vessels that supply blood to the cavernous bodies of the penis, resulting in a decreased ability to develop and maintain an erection.

Shockwave Therapy for ED of this type can be a highly effective treatment. The shockwaves are focused onto the area to be treated creating new blood vessels in penile tissue, enabling patients to achieve and maintain firm spontaneous erections.

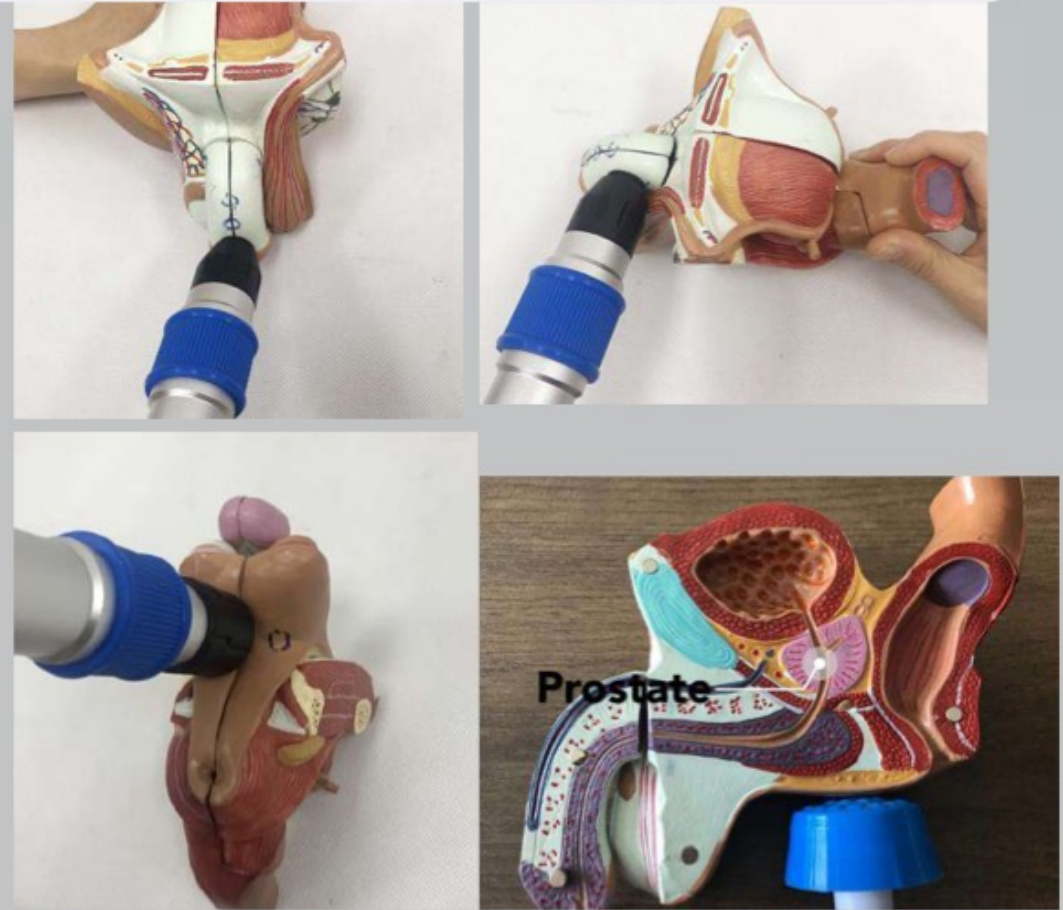

Recommended Treatment

1. Target to spongy tissue in 5 different areas
2. 300 impacts per area , total 1500 impacts in one treatment
3. Twice a week, treat consistently for 3 weeks, then stop for 3 weeks, And start next course .
4. Should be more impacts at the top end of penis, and Less impacts at the bottom of penis. Vascular-related ED problems rather than just treating the symptoms.

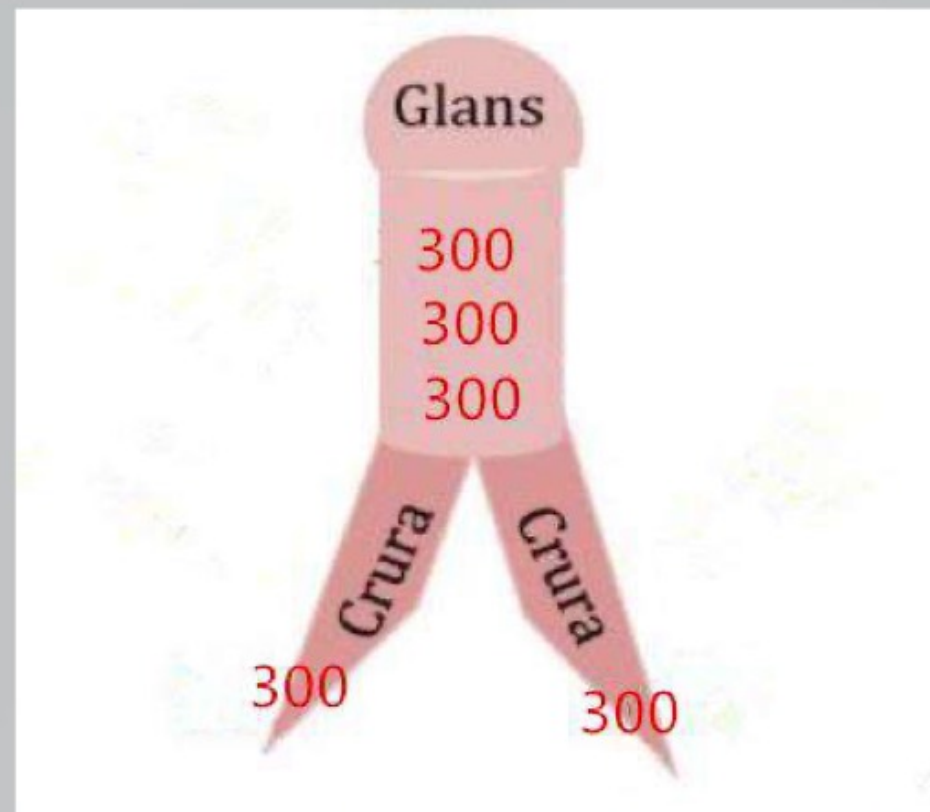

Treated Area



More Pictures




QMmed
MEDICAL AESTHETIC SOLUTIONS

Beijing QM Medical Co., Ltd.
www.zenlase.com

The Package



Our Service

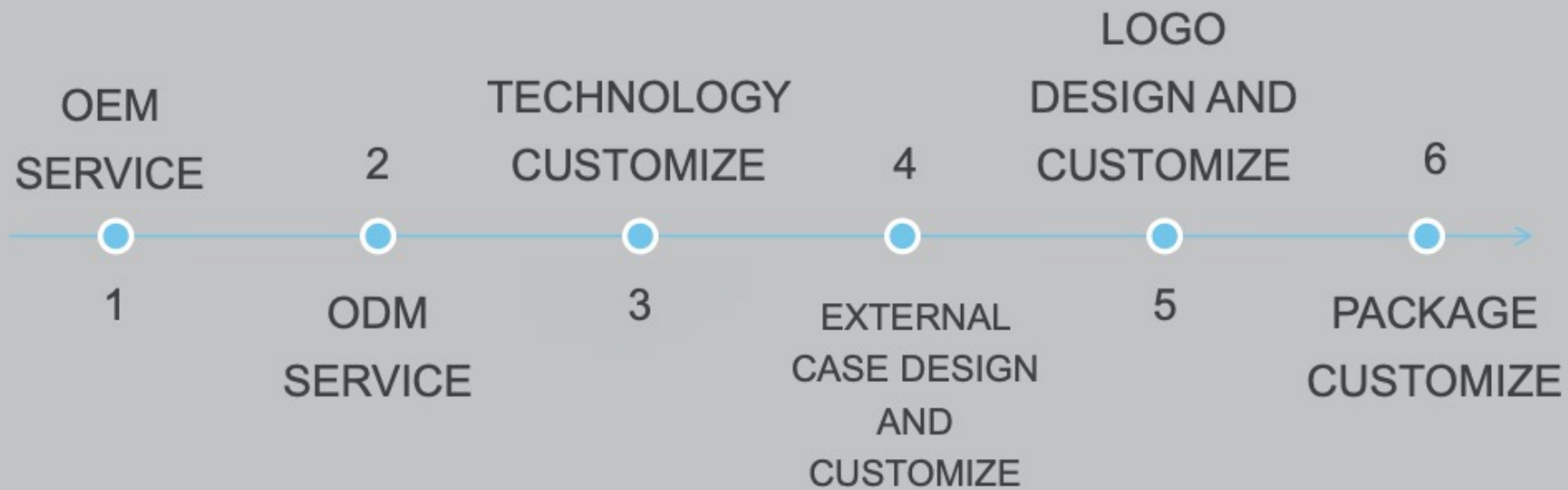

ZENLASE

Aesthetic Laser Equipment · Factory Direct

WHATSAPP

+86 151 0211 7919

EMAIL

info@zenlase.com

WEB

www.zenlase.com

HOURS

Mon-Sat 09:00-18:00 (GMT+8)

Quotations, orders, spare parts and after-sales for overseas buyers are handled by the ZENLASE team.