




INJECT NEW VITALITY INTO THE SKIN



Brighten
skin tone

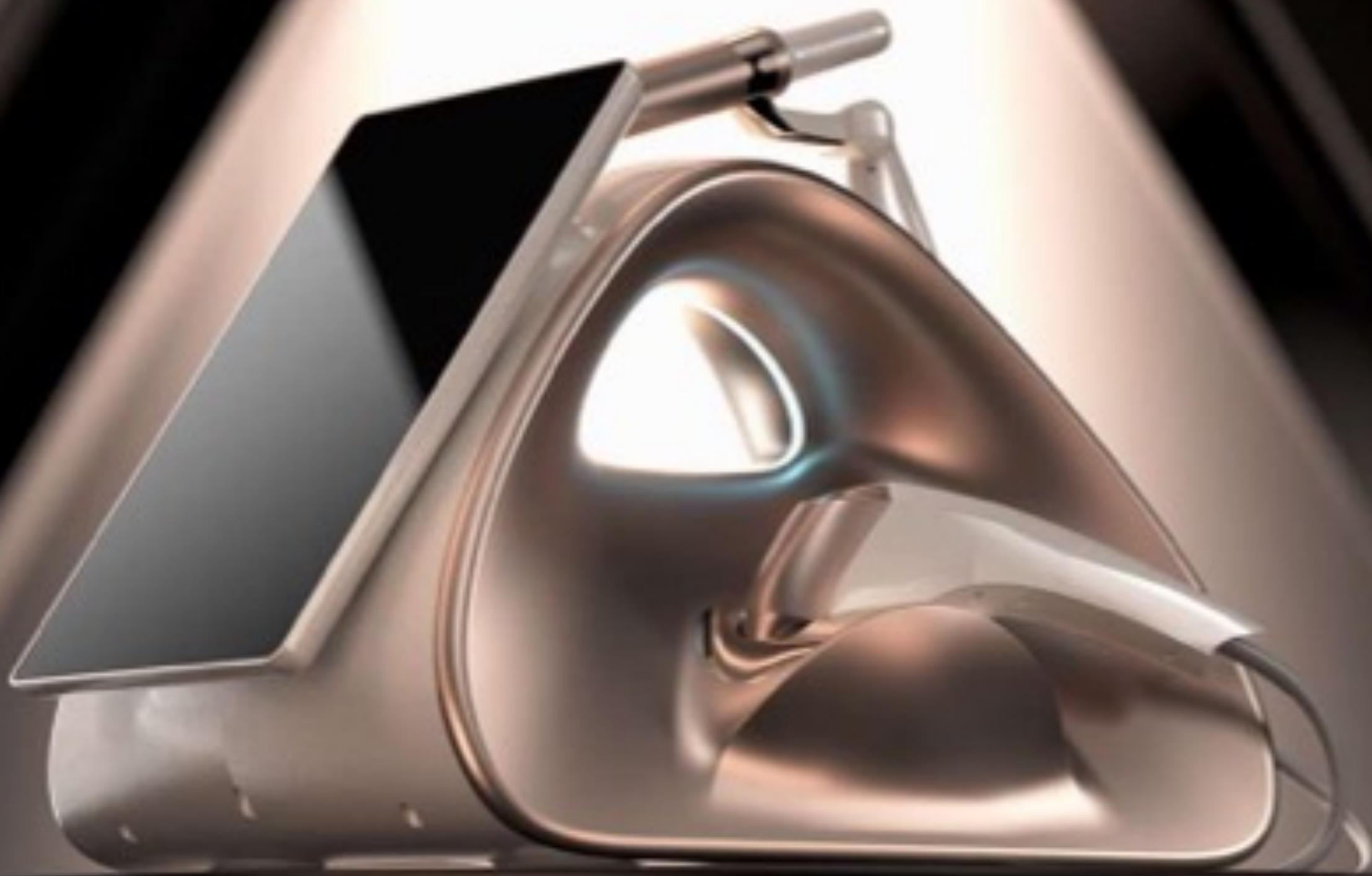

Delicate skin



Regenerated
collagen



Smooth
fine lines



Enhance
elasticity



Fat pad
reduction



Resist
gravity




Fascia
lifting

ULTRAMICRO PULSE ULTRASONIC SYSTEM

Using the ultra-micro pulsed ultrasonic system, the two-way instantaneous energy switching between the transverse linear acoustic pulse and the annular acoustic wave is realized. The single machine realizes the combination of point, line, and surface technology mode to achieve a full mode of 360 ° perfect seamless full-layer anti-aging multi-layer lifting and tightening, with precise, comfortable, fast, and durable effects.

In MP mode, the microsecond energy transmission mode has smaller and denser focus points than other instruments, higher security, comfortable experience, shorter recovery period, more uniform and more efficient point distribution.

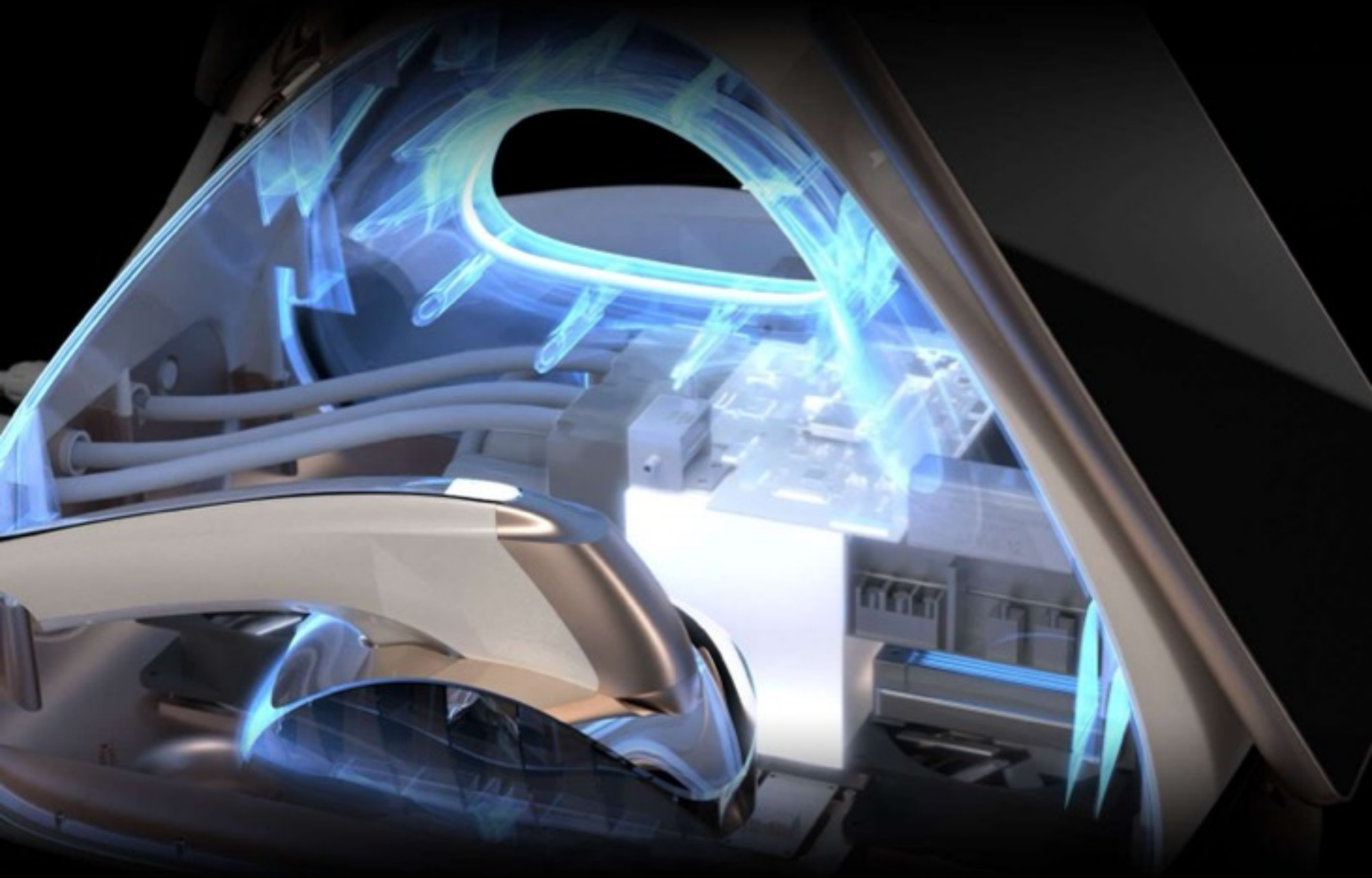

TT DTIRA

Multi-dimensional anti-aging effect

Like the “pyramid”

SAFE | STABLE | CONTINUOUS

13.3 inch

HD smart screen


Full dynamic design

human-computer interaction

Removable plug

Full sense of technology

Patented original design

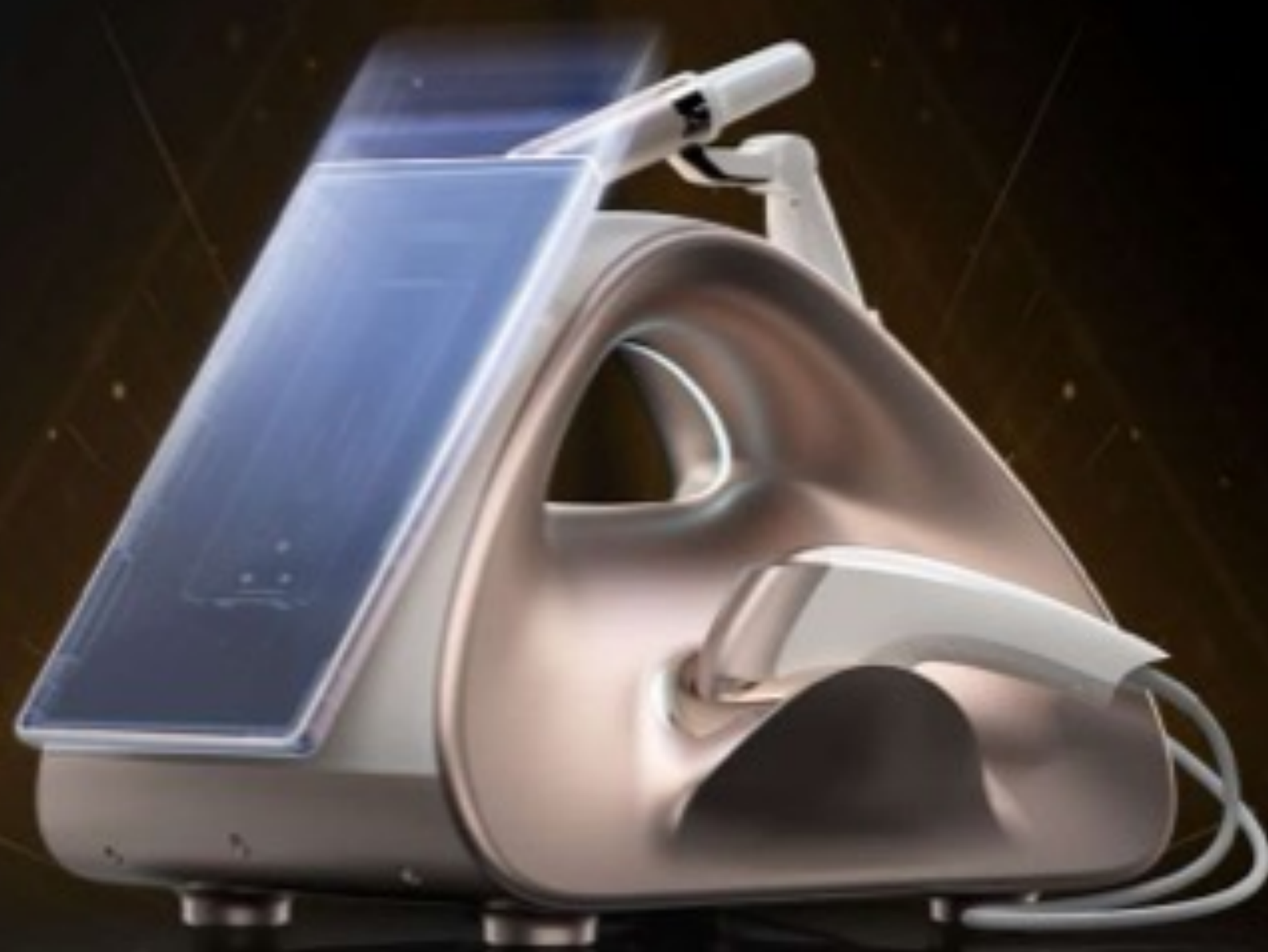



MICRO-PULSE OPERATION TECHNOLOGY&PROBE

- MP pattern (Multi-point line) [417 points/1 engine in total]
- MP pattern (Punctate) [17 points/1 engine in total]
- Mechanical effect + thermal effect Mechanical effect + thermal effect
- Two-way output innovation technology
- Multiple probes with different depths of care
- Multi-mode combined nursing

The advantages of multi-dimensional flat handle are targeted at the skin layer to the superficial fascia layer, fat layer, muscle layer, and a comprehensive customization scheme based on different skin types, operating regions, and tissue levels.

CIRCULAR OPERATION OPTIMIZATION HAND TOOL&PROBE


- MP pattern (Multi-point line) [417 points/1 engine in total]
- MP pattern (Punctate) [17 points/1 engine in total]
- Innovative technology of micro transducer
- Endless operation of the ring
- Multiple probes with different depths of care
- Multi-mode combined nursing

The world's first annular "micro-pulse ultrasound" collaborative technology, ultrasonic linear lattice focusing, can maximize the operation effect; In addition to lifting, it also brings extra firmness; Compared with traditional ultrasound mode, multiple modes can effectively reduce pain



MECHANICAL EFFECT+ THERMAL EFFECT+CAVITATION EFFECT

MP technology can more accurately transmit energy to the target skin tissue layer, generating thermal effects, tightening the fascia, activating the regeneration of collagen and elastic fibers, and achieving anti aging effects such as tightening, promoting collagen regeneration, and improving skin texture without damaging surrounding tissues.

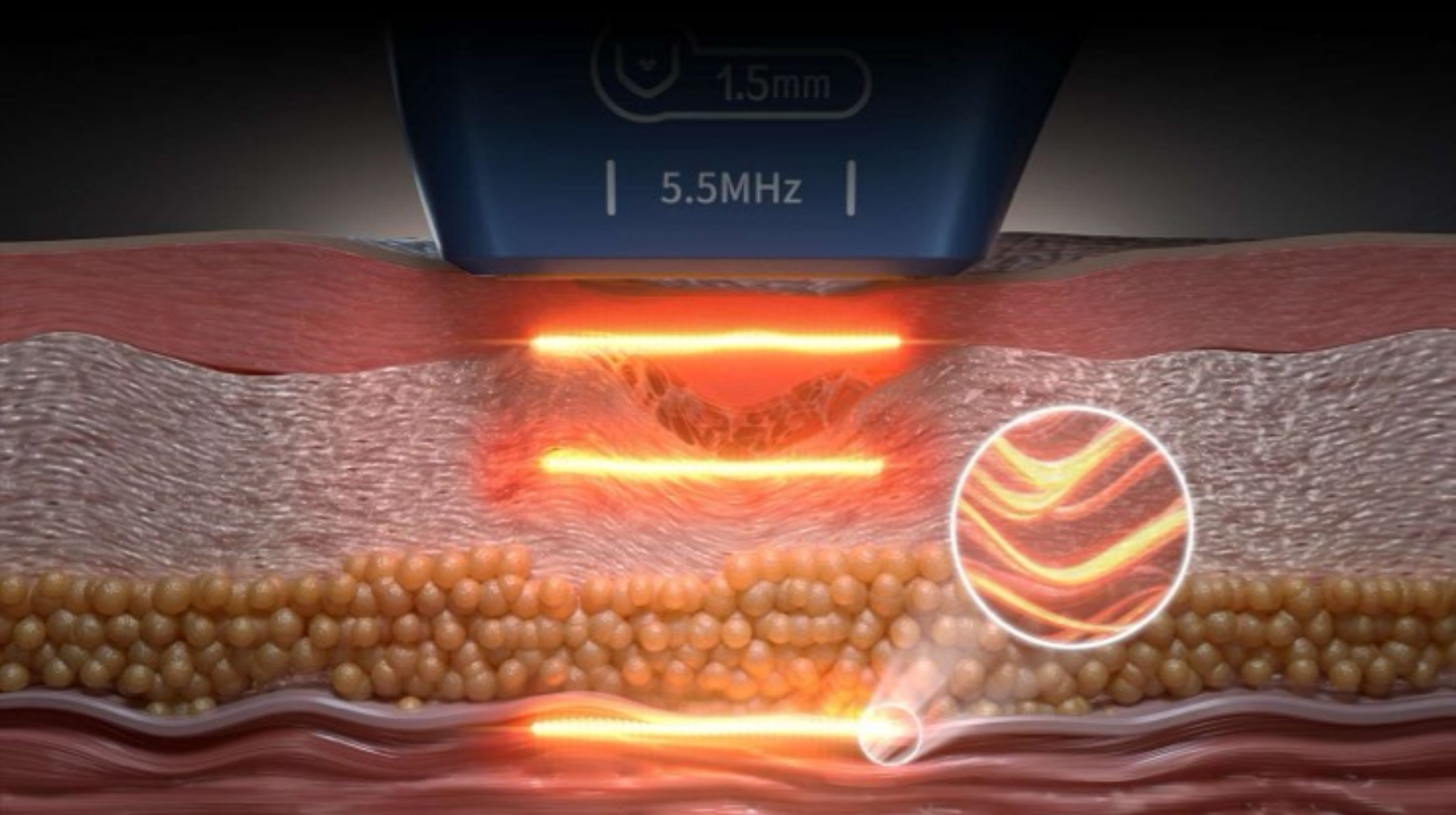

TCP IN MICRO-PULSE MODE

TCP micro-pulse mode under the micro-pulse mode realizes tiny and accurate thermal condensation point superposition, uniformly and intensively acts on all areas of the target level, and the dense and scattered energy plays a more extensive role, highlighting the overall advantages.

CIRCULAR OPERATION OPTIMIZATION HAND TOOL

The annular operating hand tool integrates the world's highest new technology ultrasonic energy output mode.

- Exquisite probe, exquisite care.
- Fusion point mode&MP mode - endless output of ultrasonic energy in circular operation.

ANNULAR RUNNING PROBE

Basic nursing effect, add new operation scheme, consolidate effect and expand adaptability Should be within the scope.



SPECIAL EYE PROBE

- 5.5 Mhz 2.0 mm special eye probe
- Unique ultra-thin probe technology
- World professional R&D application

Acting as an eye spot probe, the eye spot probe refines and solves problems such as shallow fine lines, deep genuine lines, dark circles, and eye bags.

FACE PROBE ACTS ON SKIN

Nomal (Punctate)

Micro Pulse (Linear)

- Apply gel to the skin before operation
- Operating handle vertical skin
- Do not leave the skin during operation

T - 1.5 (face) 5.5MHz


For the skin layer, superficial lines, skin color, skin texture, and pore shrinkage.

T - 3.0 (face) 5.5MHz

It directly reaches the dermis and stimulates the continuous proliferation of collagen.

T - 4.5 (face) 5.5MHz

Reach directly to the SMAS fascia layer, and tighten the muscle base to achieve tight lifting.

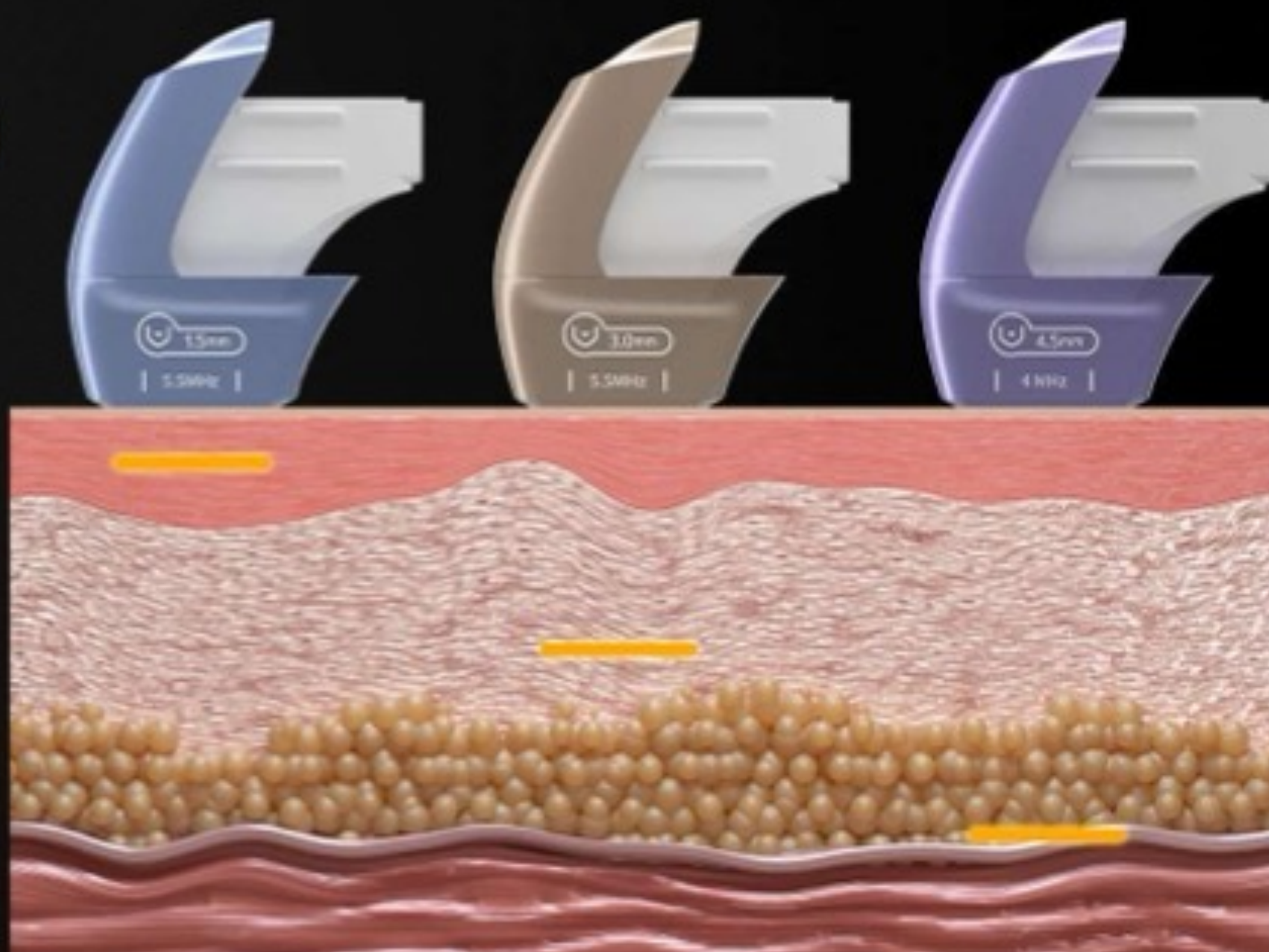

THE BODY PROBE IS APPLIED TO THE SKIN

Nomal (Punctate)

Micro Pulse (Linear)

- Apply gel to the skin before operation
- Operating handle vertical skin
- Do not leave the skin during operation

T - 6.0 (body) 2 MHz


For the skin layer, skin color, skin texture and pore shrinkage.

T - 9.0 (body) 2MHz

It directly reaches the dermis and stimulates the continuous proliferation of collagen.

T - 13 (body) 2MHz

Reach directly to the SMAS fascia layer, and tighten the muscle base to achieve tight lifting.

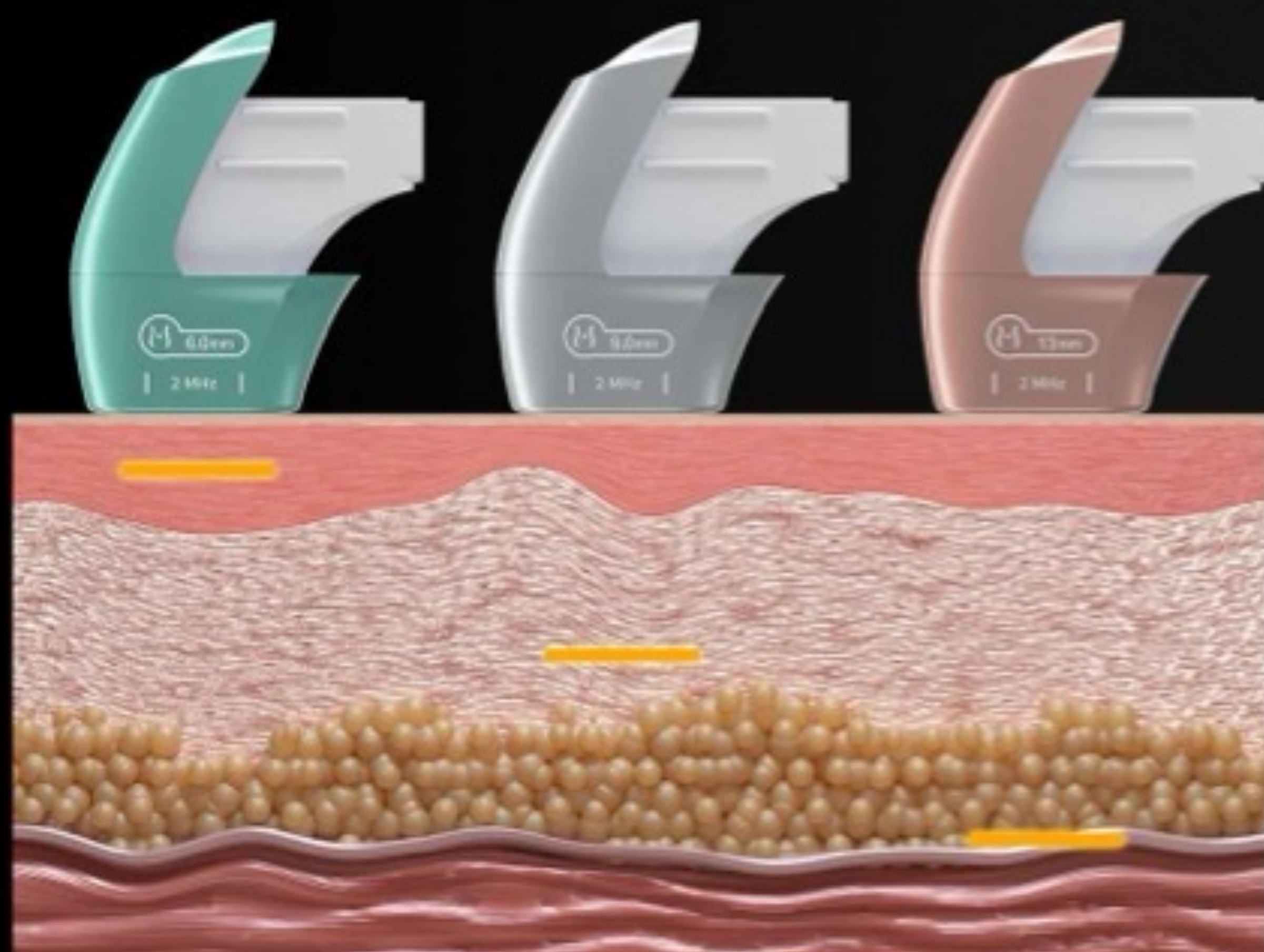

THE RING PROBE IS APPLIED TO THE SKIN



Normal
(Punctate)



Micro Pulse
(Linear)

- Apply gel to the skin before operation
- Operating handle vertical skin
- Do not leave the skin during operation

T - 1.5 (face) 5.5MHz


For the skin layer, superficial lines, skin color, skin texture, and pore shrinkage.

T - 3.0 (face) 5.5MHz

It directly reaches the dermis and stimulates the continuous proliferation of collagen.

T - 4.5 (face) 4 MHz

Reach directly to the SMAS fascia layer, tighten the bottom of the muscle to achieve a tight lift.

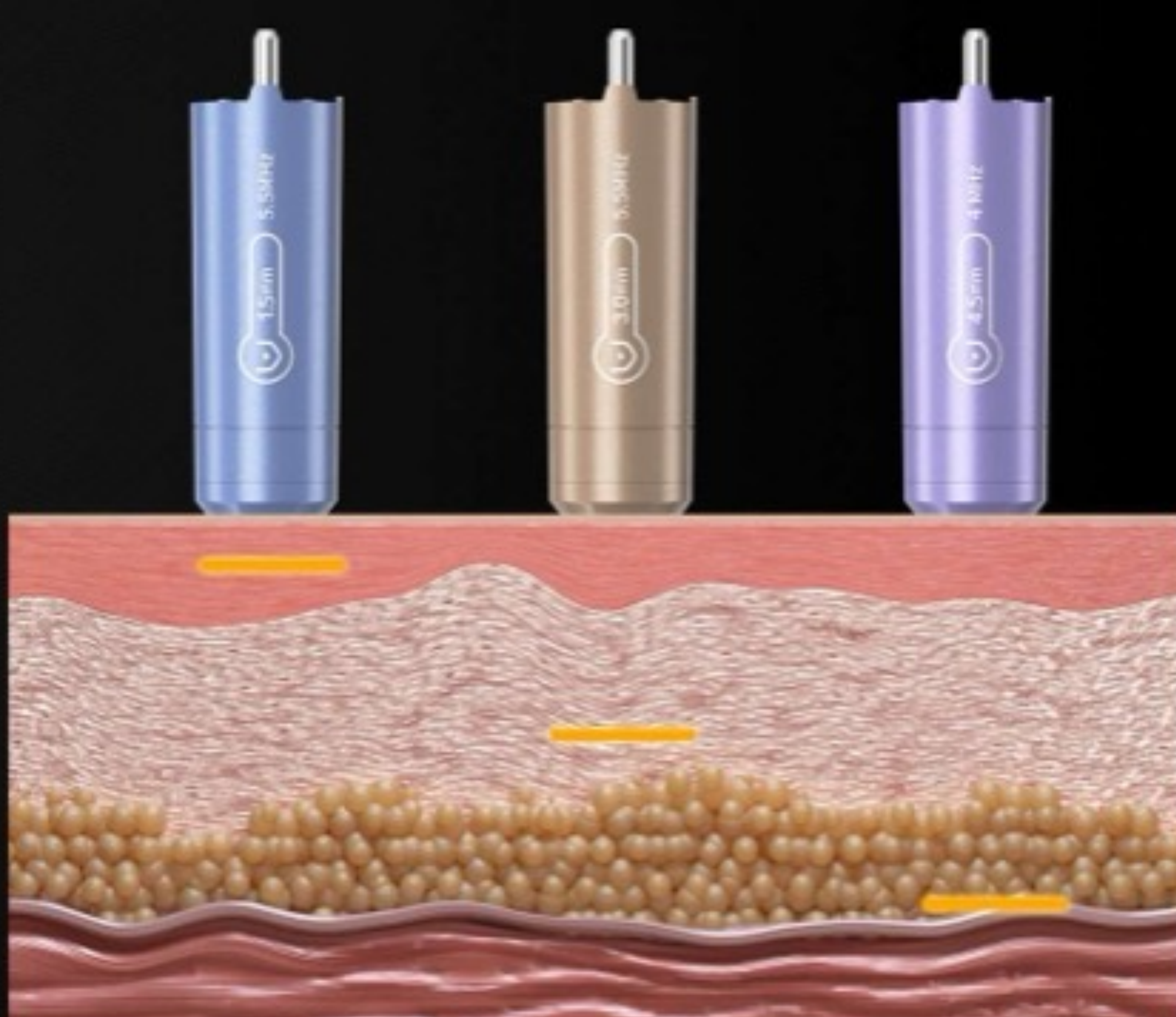


PERSONAL TAILOR FULL-LAYER ANTI-AGING

According to the age group, skin status, and customer needs, point, line, and surface can be combined to achieve a full-layer, all-dimensional anti-aging private customization program.



- Activation of epidermis and superficial dermis
- Deep dermal activation
- Instantly tighten the fascia
- Enhance the support of fascia layer and reduce the sagging tissue
- Improve cell activity and increase cell water storage capacity
- Increase facial volume, stimulate collagen regeneration, elastic fiber mesh regeneration, plump, elastic

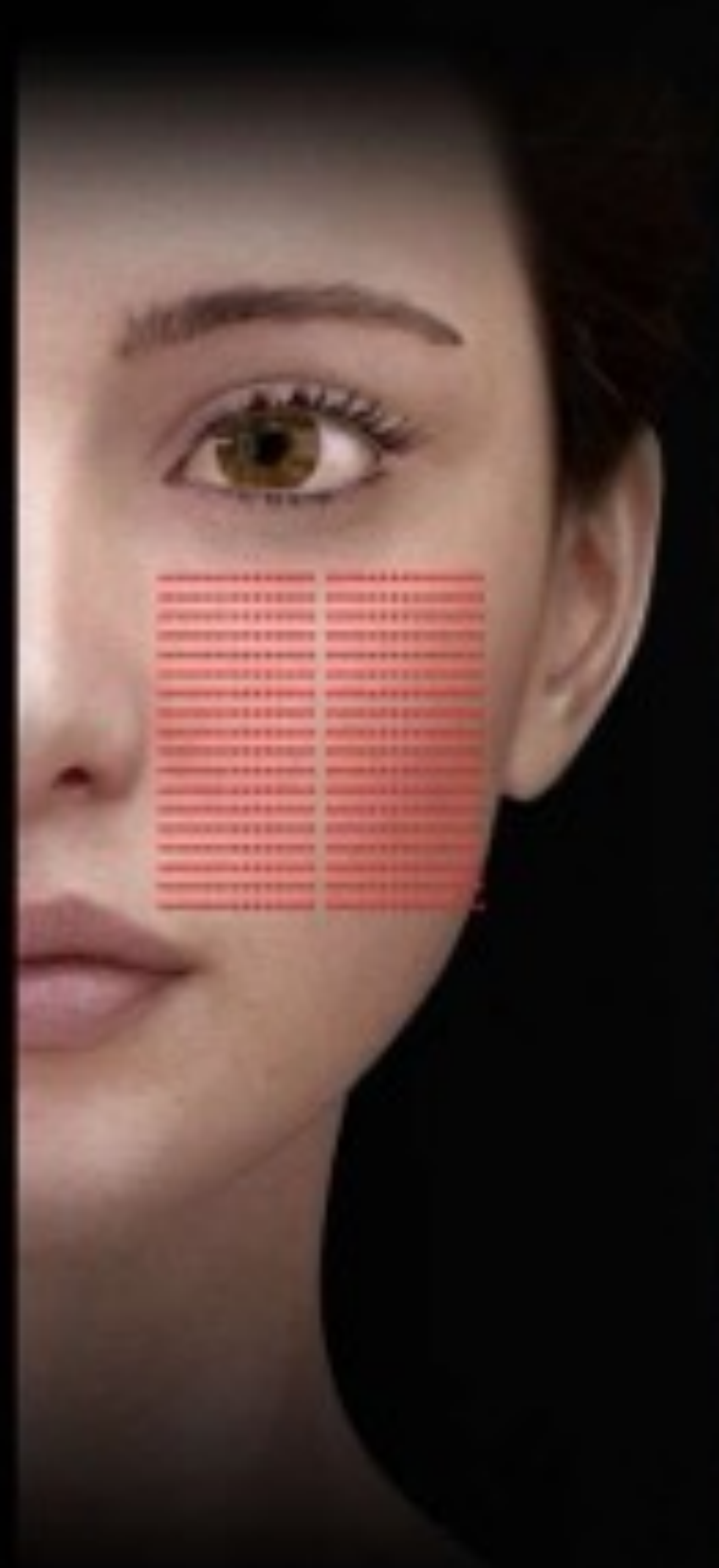

SINGLE MODE NURSING


MP mode (multi-point line formation)
[417 points/1 engine in total]

Dot mode can also be selected
[17 points/1 round in total]

The mode can be changed freely by touching the interface of the same probe

MULTI-MODE COMBINED NURSING

The micro-pulse realizes the superposition of tiny and precise thermal condensation points, and uniformly and intensively acts on all areas of the target level. The dense and scattered energy acts more widely, highlighting the overall advantages.

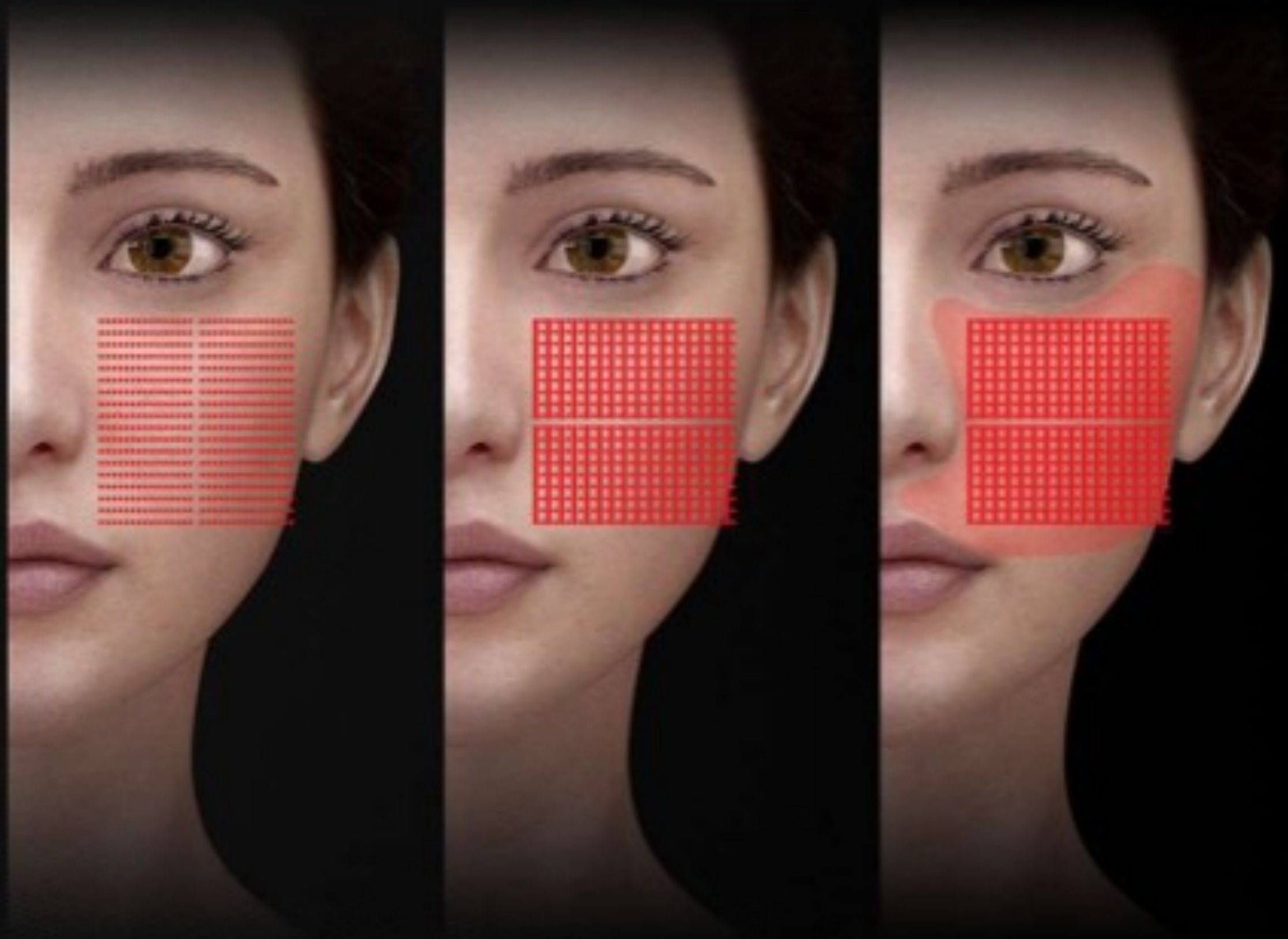


MICROPULSE OPERATION NURSING

Basic treatment efficacy, add new operation scheme, consolidate the effect and expand the scope of application

The annular MP mode can be selected as the annular point mode, and the same probe can freely change the mode by touching the interface.

PRODUCT CONFIGURATION

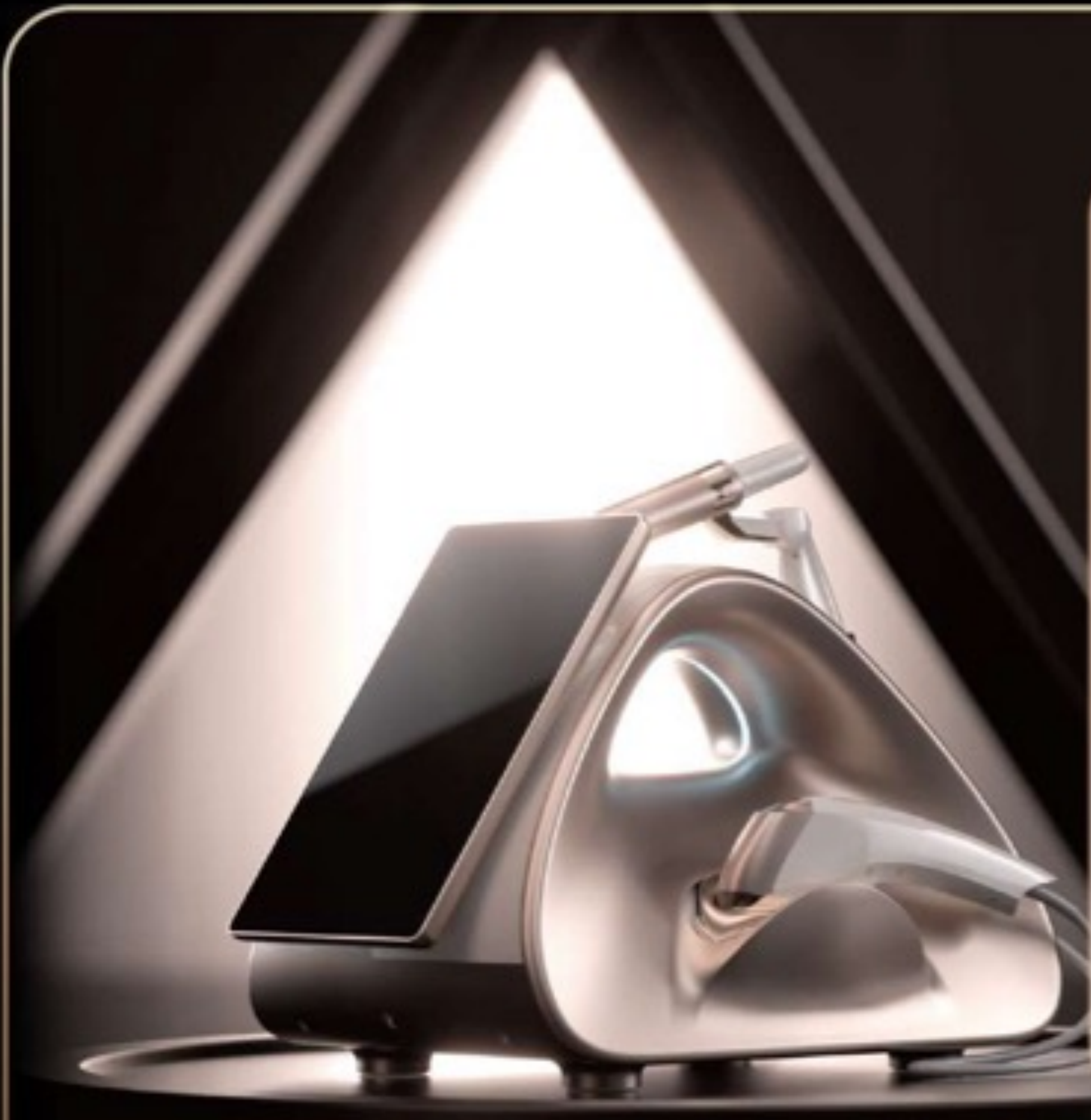

Product Name : TT MPTSTL

Model : TT MPTSTL

Power Supply : AC 110V AC 220V

Power : 650W

Frequency : 2 MHz 4 MHz 5.5 MHz

Instrument weight : 10 kg

Gross product quantity : 20 kg

Package size : 60 X 43 X 57 cm



instrument



Handle



Handle



Handle



probe



probe



Footrest



power cord



instructions



packing box

PRECAUTIONS AFTER OPERATION

- It is normal to have temporary slight redness and swelling after operation; Pay attention not to go to high temperature environment such as hot sauna yoga and strenuous exercise within 1 week after operation, and do not expose to the sun;
- Try not to make up on the day after operation (if makeup needs to be removed, pay attention to the strength, and gently remove makeup oil to clean the face);
- Do not clean your face with overheated water within a week (water that does not exceed body temperature is enough); Pay attention to water supplement. It is recommended to apply water supplement facial mask every morning and evening for 3 days or one week; Skin care and moisturizing products should be wiped 2-3 times as much as usual;
- Don't eat anything with heavy color such as soy sauce within a week, and forbid alcohol, spicy, seafood and other foods;
- It is not recommended to operate similar instruments at the same place on the same day after operation;
- After operation, it is recommended to use laser and other instruments after skin repair. (The same part shall be adjusted according to the skin condition of the day);
- There may also be a small white energy rash or skin erythema at some operating sites, which are normal, mild and temporary reactions (can be treated with ice or growth factor repair products according to the specific situation)

ZENLASE

Aesthetic Laser Equipment · Factory Direct

WHATSAPP

+86 151 0211 7919

EMAIL

info@zenlase.com

WEB

www.zenlase.com

HOURS

Mon-Sat 09:00-18:00 (GMT+8)

Quotations, orders, spare parts and after-sales for overseas buyers are handled by the ZENLASE team.